

MCKESSON

Empowering Healthcare

AcuDose – Machines

*EMT
Training*



AcuDose Machine



High Capacity



Matrix Drawer

McKesson Fridge lock



Locked Lidded Drawer



Eight full drawer slots

1

2

3

4

5

6

7

8

Which cabinet?

- You are to only use the AcuDose cabinet located between the triage rooms in the Unity Hospital Emergency Center (ED2)

User ID/password

- Your user ID will be your EMT #
- The initial password will be 1,2,3,4,5,6
- You will be prompted to change your password on initial log-in
- You will be immediately registering your biometrics
- **Do Not** share your user ID/password as this is grounds for immediate termination of your AcuDose rights

Log-In Screen

Type in your
User Name
(EMT #) &
password

AcuDose-Rx™

McKesson Automation

Version: 6.7.0

Patent Pending



User ID

PIN



Copyright © 1998 McKesson Corporation and/or one of its subsidiaries. All Rights Reserved. Protected by U.S. Patents 6,671,579 ; 6,650,964 ; 6,785,589 ; Patents Pending.

MCKESSON
Empowering Healthcare

Touch Login

Help

Drug
Reference

Training

Emergency
Dispense


Login

2D - Cardiac Care Unit


Receiving Orders

Narc Discrepancies: 1

Fast Entry – BioMetric

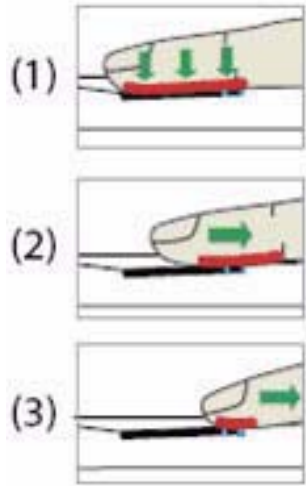

- Keep your finger flat when placing it on the black strip
- In one motion pull finger towards you
- Register four fingers (right and left/index and middle)

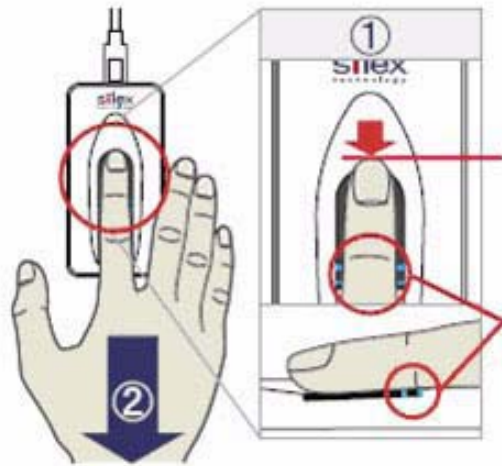

Working With Tabs

The screenshot shows a software interface for medication management. At the top, there are five tabs: **Dispense** (highlighted in green), **Return**, **Med Mgt**, **Discrepancy**, and **Admin**. A large black arrow points from the **Dispense** tab to the **Return** tab, and another arrow points from the **Admin** tab to the **Discrepancy** tab. Below the tabs is the text "Select a Patient for Medication Dispense".

The main area contains a list of patients, sorted by name. The list is as follows:

Patients sorted by: Patient Name			
ALABAMA, MASON 123458897 Unit: 4A Med-Surg Room: 410 Bed: 2			
ARIZONA, PHOENIX 123456789 Unit: 4A Med-Surg Room: 403 Bed: 1			
FLORIDA, OCALA 123456780 Unit: 4A Med-Surg Room: 404 Bed: 1			
GEORGIA, DALTON 123456776 Unit: 4A Med-Surg Room: 402 Bed: 2 Allergies!			
ILLINOIS, HINES 123456657 Unit: 4A Med-Surg Room: 401 Bed: 1			
NEVADA, VEGAS 123456784 Unit: 4A Med-Surg Room: 401 Bed: 2			
OHIO, COLUMBUS 123456875 Unit: 4A Med-Surg Room: 410 Bed: 1			

To the right of the patient list is a vertical scroll bar with up and down arrows. Below the scroll bar are several buttons: **Admit**, **Patient Info**, **Sort**, and **Logout**. At the bottom of the screen, there are three status boxes: **Training**, **Receiving Orders**, and **Discrepancies: 0**.

By touching each one of the tabs at the top it will bring up a different screen.

(EMT's will mainly need to just work with the Dispense tab.)

This could also be considered your main screen once you log in.

Selecting Patient at the Cabinet

The screenshot shows a software interface for medication dispensing. At the top, there are five tabs: 'Dispense' (highlighted in green), 'Return', 'Med Mgt', 'Discrepancy', and 'Admin'. Below the tabs is the heading 'Select a Patient for Medication Dispense'. A table lists patients with columns for patient name, unit, room, and bed. A red arrow points to the 'Dispense' tab. Another red arrow points to the 'ALABAMA, MASON' patient entry. A third red arrow points to the 'Admit' button on the right. A green box highlights a vertical scroll bar between the patient list and the buttons. At the bottom, there are three buttons: 'Training', 'Receiving Orders', and 'Discrepancies: 0'.

Patients sorted by: Patient Name			
ALABAMA, MASON 123458897 Unit: 4A Med-Surg Room: 410 Bed: 2			
ARIZONA, PHOENIX 123456789 Unit: 4A Med-Surg Room: 403 Bed: 1			
FLORIDA, OCALA 123456780 Unit: 4A Med-Surg Room: 404 Bed: 1			
GEORGIA, DALTON 123456776 Unit: 4A Med-Surg Room: 402 Bed: 2 <i>Allergies!</i>			
ILLINOIS, HINES 123456657 Unit: 4A Med-Surg Room: 401 Bed: 1			
NEVADA, VEGAS 123456784 Unit: 4A Med-Surg Room: 401 Bed: 2			
OHIO, COLUMBUS 123456875 Unit: 4A Med-Surg Room: 410 Bed: 1			

Under the Dispense Tab at the top, select the patient.

If your patient's name does not appear in the list:

Hitting the admit button will allow you to admit someone at the station.

ADMIT SCREEN

LOPEZ, OLIVIA Details

Admit | Optional | Misc | Allergies

Patient Information

Site Patient ID: 123456783

Name: LOPEZ, OLIVIA

Connect-Rx Patient ID: 000000002356

Admit Date: 03/23/2007 07:00 Class: Inpatient

Del. Site: Primary delivery site Source: External

Bed: 4A 402 1 Station: INT1

Status: Admitted Verified: Yes

Cancel

5C - MedSurg Receiving Orders **Narc Discrepancies: 4**

It will prompt you to Enter: Name (last, first) Financial number (if available)


If code 500, add 500 and today's date for the name

Selecting the Medication from the list

Dispensing for Patient: LARSEN, DAWN

Brand Generic Dispense

Medications sorted by: Generic Name

100 GRAND 20 MCG BAR HEALTH <i>In Stock: 0</i>	Out of Stock
ACETAMINOPHEN 1 TAB CAPSULE ACETAMINOPHEN 500MG CAPLET <i>In Stock: 0</i>	Out of Stock
ACETAMINOPHEN-COD 30 MG TABLET APAP-CODEINE #3(ALL SCRIPS) <i>In Stock: Count Required</i>	
 ALBUTEROL 17 GM INHALATION VENTOLIN <i>In Stock: 8</i>	
ALBUTEROL 3 ML SOLUTION ALBUTEROL SOLN PREMIX 0.083 <i>In Stock: 3</i>	
ALPRAZOLAM 0.5 MG TABLET XANAX <i>In Stock: Count Required</i>	
ALUMINA MAG AND SIMETHICONE 30 ML LIQUID MAG AL PLUS UD <i>In Stock: 10</i>	
AMOXICILLIN-CLAV 500 MG TABLET AUGMENTIN 500MG <i>In Stock: 10</i>	

Kit Detail
Drug Reference
Cancel

LEARNING CENTER Receiving Orders **Narc Discrepancies: 1**

Having the drug highlighted in Green indicates that is the med you have selected to take out. (You can take out more than one drug at a time, by selecting several drugs at once.) If you decided you do not want that drug you can touch it and it will be unselected.

*Green highlight means the item is selected.

Dispense Screen.... Continued

By hitting this button it will allow you to dispense the highlighted Green Items.

These arrows allow you to move the list of drugs up and down. You can also start typing the name of the drug and it will bring up a text search bar above the list of meds.


Backs you out to main screen.

Dispensing for Patient: LARSEN, DAWN

Brand Generic

Dispense

Medications sorted by: Generic Name

100 GRAND 20 MCG BAR HEALTH <i>In Stock: 0</i>	Out of Stock
ACETAMINOPHEN 1 TAB CAPSULE ACETAMINOPHEN 500MG CAPLET <i>In Stock: 0</i>	Out of Stock
ACETAMINOPHEN-COD 30 MG TABLET APAP-CODEINE #3(ALL SCRIPS) <i>In Stock: Count Required</i>	
 ALBUTEROL 17 GM INHALATION VENTOLIN <i>In Stock: 8</i>	
ALBUTEROL 3 ML SOLUTION ALBUTEROL SOLN PREMIX 0.083 <i>In Stock: 3</i>	
ALPRAZOLAM 0.5 MG TABLET XANAX <i>In Stock: Count Required</i>	
ALUMINA MAG AND SIMETHICONE 30 ML LIQUID MAG AL PLUS UD <i>In Stock: 10</i>	
AMOXICILLIN-CLAV 500 MG TABLET AUGMENTIN 500MG <i>In Stock: 10</i>	

Kit Detail

Drug Reference

Cancel

LEARNING CENTER Receiving Orders Narc Discrepancies: 1

Dispense Screen

The screenshot shows a software interface for entering a quantity to dispense. At the top is a green header with the text "Enter Quantity to Dispense". Below this is a large input field labeled "Dispense Quantity" containing the number "0". To the right of this field are plus and minus buttons. A numeric keypad is located in the center, with buttons for digits 1-9, 0, and a "Backspace" button. To the right of the keypad is an "Open Drawer" button. At the bottom of the screen, there is a section for drug information, including "ACETAMINOPHEN 1 TAB CAPSULE", "ACETAMINOPHEN 500MG CAPLET", and "LOPEZ, OLIVIA". Below this, it says "Order: ACETAMINOPHEN 1 TAB CAPSULE" and "Unable to determine the last dispense." At the very bottom, there are three status bars: "2D - Cardiac Care Unit", "Receiving Orders", and "Narc Discrepancies: 1".

Annotations with red arrows point to the "Dispense Quantity" field, the "Open Drawer" button, and the drug information section.

This is where you will type the quantity that you would like to dispense.

It will then prompt this "Open Drawer button" to light up

The bottom will show the drug you are dispensing.

How to Clear Cabinet Error

- Either:
 - Log out of the cabinet and log back in
 - Go to “Admin Tab” – Hit Clear Error Button
- Contact EC charge nurse or pharmacy with issues

Congratulations you have completed the EMT Training!!

- You now need to review the Medication Administration Agreement

By reviewing this PowerPoint and signing the agreement you have a full understanding of the proper way to use the Acu-Dose Machine at Unity Hospital.

If you do not understand please seek help from EC Charge Nurse or Pharmacy.

THANK YOU!